

PATENTTI- JA REKISTERIHALLITUS  
NATIONAL BOARD OF PATENTS AND REGISTRATION

PCT/FI 03 / 00247

Helsinki 20.5.2003

ETUOIKEUSTODISTUS  
PRIORITY DOCUMENT

REC'D 10 JUN 2003

WIPO

PCT

Hakija  
ApplicantFocus Inhalation Oy  
TurkuPatenttihakemus nro  
Patent application no

20020607

Tekemispäivä  
Filing date

28.03.2002

Kansainvälinen luokka  
International class

A61K

Keksinnön nimitys  
Title of invention"A method for treating carrier particles and its use"  
(Kantoinen hionta)

Täten todistetaan, että oheiset asiakirjat ovat tarkkoja jäljennöksiä  
Patentti- ja rekisterihallitukselle alkuaan annetuista selityksestä,  
patenttivaatimuksista ja piirustuksista.

This is to certify that the annexed documents are true copies of the  
description, claims and drawings originally filed with the Finnish Patent  
Office

**PRIORITY DOCUMENT**  
SUBMITTED OR TRANSMITTED IN  
COMPLIANCE WITH  
RULE 17.1(a) OR (b)

Marketta Tehikoski  
ApulaistarkastajaMaksu 50 €  
Fee 50 EUR

Maksu perustuu kauppa- ja teollisuusministeriön antamaan asetukseen 1027/2001 Patentti- ja rekisterihallituksen maksullisista suoritteista muutoksineen.

The fee is based on the Decree with amendments of the Ministry of Trade and Industry No. 1027/2001 concerning the chargeable services of the National Board of Patents and Registration of Finland.

Osoite: Arkadiankatu 6 A  
P.O.Box 1160  
FIN-00101 Helsinki, FINLANDPuhelin: 09 6939 500  
Telephone: + 358 9 6939 500Telefax: 09 6939 5328  
Telefax: + 358 9 6939 5328

A method for treating carrier particles and its use.

The present invention relates to a method for treating a particulate carrier for an inhalation powder improving stability and flow properties of the carrier. The invention further concerns the carrier and a pharmaceutical preparation for inhalation purposes containing said carrier.

Micrometer-millimetre size particulate materials, organic or inorganic, are normally not spherical nor rounded but rather edged or rough surfaced after crystallisation or other manufacturing methods.

The particle characteristics of the raw materials strongly affect the final quality of the pharmaceutical product and therefore very strict requirements are applied on these properties in the pharmaceutical industry. A special case of such material science in pharmaceutical industry is the formulation for a powder inhaler. The overall performance of the powder inhaler is highly dependent on the characteristics of the powder components, e.g. particle size distribution, crystal morphology, shape and surface roughness of the particles and interparticle attraction forces, including static charges. An important aspect which must be taken into account with handling and processing of inhalation powders is the prevention of microbial contamination.

Delivery of drugs into the lungs is common in asthmatics and is usually done with a metered dose inhaler (MDI; pressurised aerosol) or a powder inhaler. Irrespective the dosing device, the delivered particles must be no larger than 5 micrometers in respect of the tested MMAD (mass median aerodynamic diameter) in order to be able to deposit in the deep lung. It is expected that the use of powder inhalers will strongly increase and systemic delivery of new drugs, including large molecules, will be a marked target for powder formulation technology.

## 2

The classical powder formulation consists of an air-jet milled micronsize drug and a carrier sugar, commonly lactose monohydrate. The mean particle size of the carrier is mostly between 50 and 100 microns and the particle size distribution is broad resembling a Gauss' curve. The maximum size is about 300 microns, larger may cause irritation in upper airways. The form of commercial lactose is typically an elongated triangle, called "tomahawks". Glucose can be used as carrier and the particles do not differ markedly from lactose in respect to the particle size distribution or particle form. Mannitol, sucrose and trehalose are under investigations for carrier sugars. Classified samples are offered for test purposes by some sugar manufacturers and such carrier materials will be available for innovators in manufacturing scale.

The drug content in a powder formulation is typically less than 10 w-%. It is calculated that with a drug concentration of 5-10 % the carrier particles are covered by a uniform layer of drug particles. Higher concentrations may cause segregation and will demolish the flow properties of the formulation. Excellent flow properties are extremely important for multiple dose powder inhalers (MDPIs), where each dose must be accurately metered by manoeuvres done by the patient. Capsule filling of unit dose devices or filling of the blisters for blister-based devices can be done with formulations having poorer flow properties.

There are some main factors which regulate the performance of the formulation in use:

### 1. Force of attraction between drug and carrier particles

Micronsize particles tend to adhere very strongly to each other and carrier particles. During inhalation the drug particles should be liberated again in order to be able to go into the deep lung. This disintegration takes place due to turbulence, shear and centrifugal forces. If the carrier particles are rough

and edged, drug particles may be hidden in the high-energy corners and cavities and they are resistant to shear and turbulent flow. The degree of loose drug particles can be easily metered with a cascade impactor, where the fine particle fraction (FPF %) and mass median aerodynamic diameter (MMAD) of the delivered drug particles can be calculated.

Van der Waals forces are important in respect to the adherence when small distances are concerned. If two particles have a direct contact with large contact area, the attraction is high and they are difficult to separate. If there is a primary layer of very small particles on the carrier surface, the secondary drug particles tend to adhere relatively loosely on the carrier surface due to diminished Van der Waals forces. It means that one can increase the FPF % if very small particles are mixed to the formulation, preferably before the drug particles. Ideally, the material is the same which is used as the coarser carrier. In literature, the particle size of added small particles has been about the same as the particle size of micronised drug.

If the surface of any of the formulation components is modified, the performance of the product will change due to altered particle-particle interactions. This refers also to particle mean size and particle size distribution of the components.

## 2. Physical stability of the components

The components should be physically stable or in their thermodynamically lowest energy level. If not so, the component will change its physical state more or less slowly, accelerated by increase in temperature and humidity. The change is seen as altered performance and is a common reason for impaired shelf-life of the product. Air-jet milling of the drug creates easily amorphous material on the surface of drug particles. Vigorous dry mixing may do the same to all components. The formation of amorphous material is

highly drug specific. Some drugs may transform into a totally amorphous state whereas some do not change at all. It is commonly believed that the amorphous content in the micronised drug is mainly responsive for the impaired physical stability of inhalation powders. The role of the carrier has remained more unclear in this respect.

### 3. Factors affecting dose accuracy

Concerning MDPIs, the accuracy of the dose metering mechanism of the device is decisive. In most cases the dose is metered to a dose slot or slots to be transferred to inhalation air stream. This volumetric dose metering may work accurately only, if the formulation shows proper and unaltered flow properties through the shelf-life. If the formulation is not physically stable, changes in the morphology may cause agglomeration of the powder, followed by impaired flow properties and dose accuracy. If there are too many micronsize particles (more than 10 w-%), in the formulation, the flow properties may initially be impaired and the formulation is even more sensitive to further disturbances, e.g. to unstability of some component. In respect to flow properties, the sensitivity of the dose metering system may vary between different MDPIs.

Ideally, the best dose accuracy and longest shelf life is obtained if the formulation is physically stable, it is protected against moisture with a desiccant and the flow properties remain unchanged in actual use circumstances.

### The invention

It was discovered that edged and rough carrier sugars could be rounded and polished by treating a carrier suspension some hours with a high performance disperser. By introducing the suspension through a classifying filter pack during the procedure a certain range of particle size of polished sugar

## 5

could be obtained. Drug formulation made of polished carriers showed enhanced performance and characteristics, especially in physical stability, when used in a MDPI.

- 5 There are some known methods to modify the shape of carriers. In US  
6,153,224 classified lactose particles were gently milled in a ball-mill in order  
to remove the asperities on the surface as small grains, which will be reat-  
tached on the high-energy sites of the carrier. When a small amount of a  
ternary agent (L-leucine) was milled with lactose, the grains of the agent  
10 adhered to lactose. As a result, such a modified carrier showed decreased  
adhesion to drug particles and hence, better small particle fraction in labora-  
tory tests was achieved.

- 15 Iida et al. (Chem. Pharm. Bull. 49(10) 1326 1330 (2001) Vol. 49, No. 10)  
removed protuberances from the surface of lactose carrier by controlled dis-  
solution. The resulting particles were rounder and without sharp edges com-  
pared to untreated lactose. Drug mixtures made of such lactose showed im-  
proved flow properties and better fine particle fraction, when packed in cap-  
sules and used in a powder inhaler.

- 20 No methods concerning carrier polishing with a high energy disperser or a  
corresponding mixing device based on vigorous mixing of a suspension has,  
however, been described.

- 25 The patent application WO 02/00197 A1 Staniforth et al. discloses a method  
for making microfine composite particles. This is done preferably by wet-  
milling the components in a ball mill. It is also mentioned that a high-energy  
liquid homogenizer can be used for the purpose. In this connection, polish-  
ing of larger carriers is not mentioned.

Microencapsulation of carriers by spray drying and some other methods for coating carriers have been reported, but these include no polishing or abrasion of the carrier surface.

- 5 To improve the separation of carrier and active particle the invention suggests that carrier is abraded suspended in a liquid medium into which the carrier is essentially insoluble, the liquid medium is removed and the carrier recovered. The so abraded or smoothened carrier particles have been found to more efficiently deliberate the active particles adhered to the carrier. Also  
10 the physical stability of the treated carrier is enhanced. The flow properties of the treated, filtered and dried carrier were clearly improved.

- The abrasion is preferably conducted with a mixer device such as a high performance disperser using an effect below that required for crushing the carrier particles thus avoiding breaking up the particles to be treated. Preferably  
15 the abraded carrier is at least partly covered with fine particles.

- The invention also concerns a carrier for an inhalation powder, which carrier is stable and possesses good flow properties, characterized in that the carrier  
20 is abraded suspended in a liquid medium, in which said carrier is essentially insoluble.

- A further feature of the invention is a preparation for inhalation purposes comprising an active agent, a carrier and optional excipients used in inhalation preparations. The carrier in this preparation is at least partly abraded  
25 suspended in a liquid medium, in which the carrier is essentially insoluble. An especially advantageous preparation contains in addition to the abraded carrier also a micronised carrier. Such preparation has even a more prolonged shelf-life than a preparation manufactured of drug and polished carrier alone.

## Description of the test methods.

The first experiments were done by treating the n-hexane (Mallinckrodt Baker BV, the Netherlands) suspension of Pharmatose® 325 M lactose monohydrate (DMV, The Netherlands), mean particle size 60 microns, for some hours with an Ultra-Turrax® high performance disperser IKA T 25 Basic (20.000 rpm) (IKA GMBH & Co KG) in a decanter. The batch size was some tens of grams. It was found that up to 30 % of the initial amount of lactose was abraded to micronsize particles, which could be filtered away. In the next step the disperser was provided with a flow-chamber, suspension vessel and with an ice-bath cooled recirculation line. The product was obtained by filtering the treated suspension through a 40 micron filter, followed by vacuum drying.

The pilot-scale polishing equipment is based on IKA SD 41 Super-Dispax® high performance disperser (IKA GMBH & Co KG), equipped with a flow-chamber for circulation of the feed suspension. An on-line filter block was used to separate smaller than 10 micron particles as waste and return the larger particles to the flow-chamber. The principle of the filter is explained further in US 6,027,656. The filter principle enables to recirculate or waste more than one main range of particle size, if more than two filter planes are used.

The system is illustrated schematically in figure 1 presenting a pilot scale polishing equipment.

A water cooled suspension vessel 2 is equipped with a mechanical mixer 1 and below a Super Dispax® water cooled stator/rotor chamber 3 equipped with a motor 11. The suspension obtained in vessel 2 is fed to the chamber 3 (flow-chamber) for the rotor/stator treatment and then to the filter device 4 with a motor 10 and two filters, the first one being a coarse filter 5 (pore size



40  $\mu\text{m}$ ) and the second one is a fine filter (pore size 0,5  $\mu\text{m}$ ). An abraded screened product is obtained from filter output 9. Coarse fraction 7 and the very fine fraction 8 are returned to vessel 2.

- 5 Three different liquids were tested as lactose suspension: n-hexane (Mallinckrodt Baker BV, the Netherlands), 2-propanol (Mallinckrodt Baker BV, the Netherlands) and a mixture of non-flammable perfluoroethers [Galden® (Aurimont, Italy)]. Even if the densities of the liquids were very different (n-hexane = 0.7 g/cm<sup>3</sup>, 2-propanol 0.8 g/cm<sup>3</sup>, Galden® 1.6 g/cm<sup>3</sup>), no difference in the polishing efficacy was found.

A validation program of the system was carried out with Pharmatose® 110 M lactose monohydrate (DMV, the Netherlands) and 2-propanol as the suspending liquid. The studied parameters were lactose:2-propanol ratio, roundness (shape factor), RPM of the rotor and the distance of stator/rotor. The coarse filter was 10  $\mu\text{m}$  and the fine filter 0.5  $\mu\text{m}$ . The time of the treatment was 3 hours.

The results revealed that only RPM of the rotor is of practical importance. In practice, 100g carrier/1 litre produces a good suspension. Distance of stator/rotor should be at least twice the diameter of the largest particles. If not, the particles are crushed and not polished. Results on the effect of RPM on particle size (10%, 50% and 90% with smaller size) and roundness are seen in Figure 2.

25

30



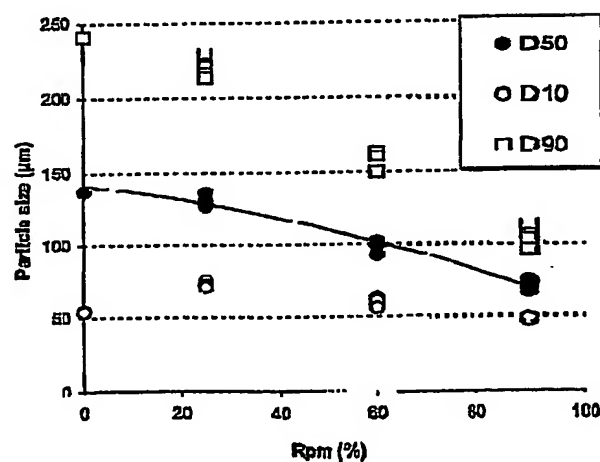


Figure 2. Mean particle size and shape factor in function of RPM in polishing 110 M lactose

- 5 The desired particle size may be obtained by choosing the appropriate rotor/stator distance and/or rotation speed of the mixer. Said distance is material specific, for greater particles a distance of 0,5 mm may be sufficient, whereas the distance may be reduced if smaller sized particles are treated. The higher the rotation speed is the smaller is the resulting average particle size. Finally, the treatment time is decisive.
- 10

- One should use at least 60 % of the maximum rotor speed (13.000 RPM) and at least 80 %, if maximum roundness is wanted. However, if also maximum polishing is wanted, one should use first RPM of at least 80 % and after that 60 hours with 60 % or less, when the final polishing is done gently and gives smoother surface. The result of a routine treatment of 110 M lactose in propanol (3 h, 80 % RPM) is seen in figure 3.
- 15

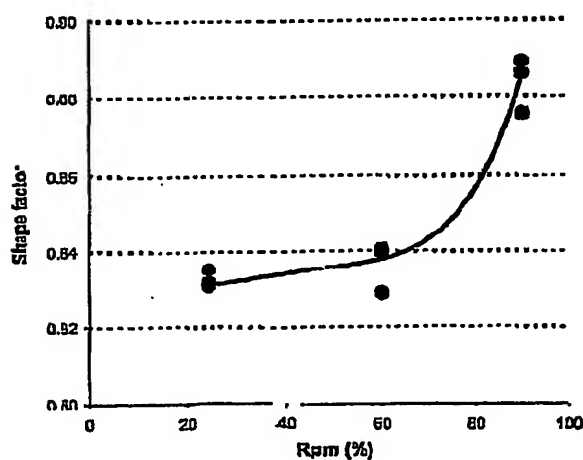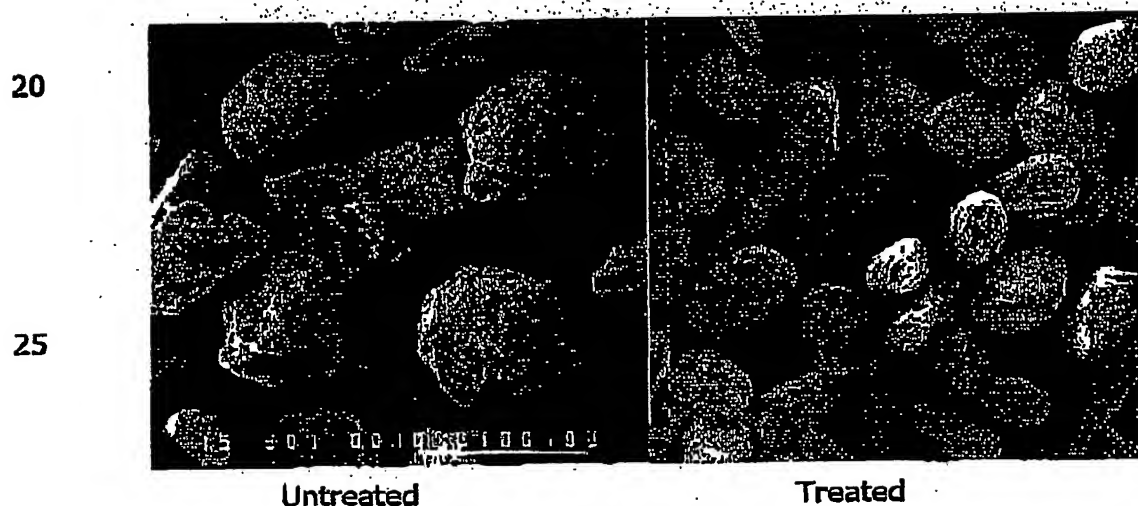

30 Figure 3. Result of a routine treatment

10

Another example is explained, when 110 M glucose anhydrous (Kirsch Pharma GmbH, Germany) was treated in Galden® 100 (Ausimont, Italy): Ultra Turrax® basic provided with a flow-chamber and ice-bath cooled recirculation was used as a disperser. Treatment time was 1.5 h, disperser speed was 22.000 RPM and the amount of glucose anhydrous was 150 g. 1500 ml of GALDEN® 100 was used as medium.

After the treatment the suspension was filtrated through a 40 micron filter and was washed several times with n-hexane in order to remove the residual small particles. Then the filtered mass was dried in vacuum. The dry product was sieved through a 150 micron sieve.

Microscope photos and particle size distributions of the starting material and the final product are seen in Figure 4.

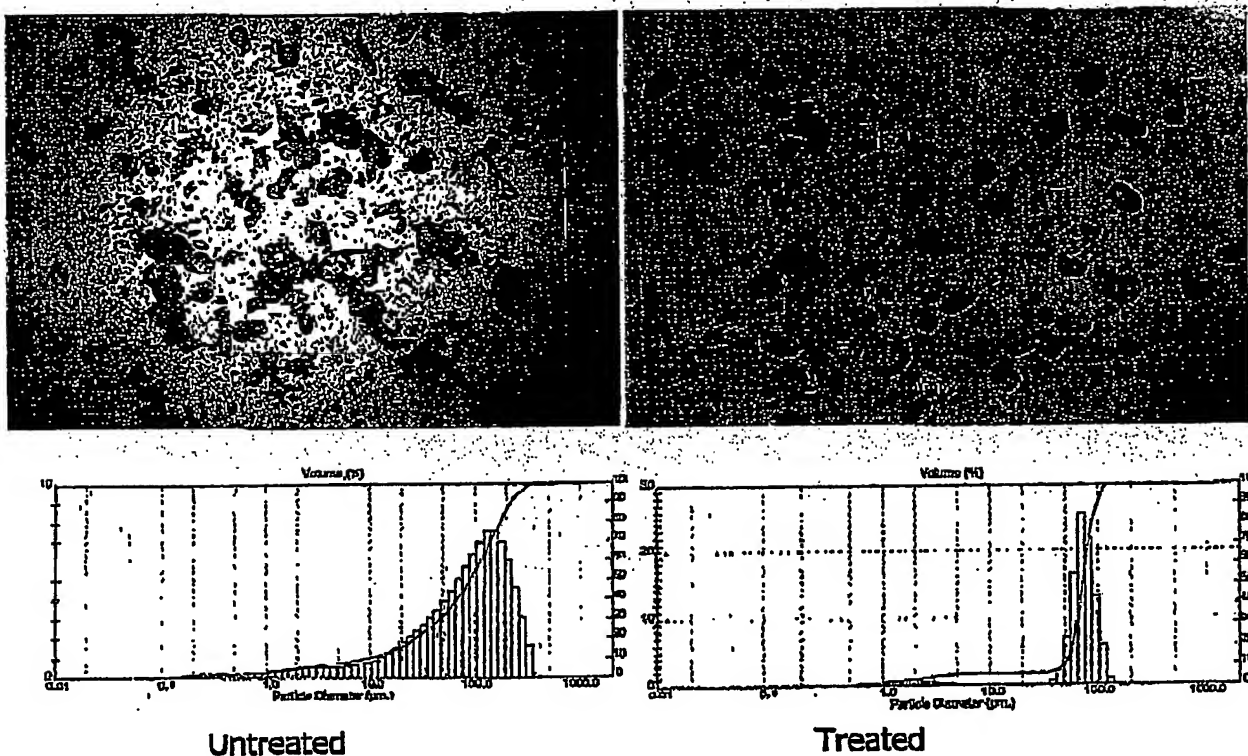


Figure 4. Microscopic image and particle size distribution of untreated and polished glucose

According to microcalorimetric studies, untreated lactose contained detect-  
5 able amounts of unstable (amorphous) material. When the study was re-  
peated with polished lactose of the same manufacturing batch, no signs of  
amorphous material were found. It is obvious that the amorphous matter  
was located at the surface of lactose and was removed by polishing. Lack of  
10 amorphous substance on the surface of polished carrier particles is most ob-  
viously the reason for enhanced stability of the final formulations. The im-  
provement in stability was surprisingly clear and indicates the importance of  
the carrier in addition to the micronised drug in respect to physical stability  
of the formulation.


#### 15 Test results with formulations containing polished carrier

The tests were carried out as follows: The formulations were manufactured  
by the wet-mixing method explained in Finnish patent No. 105.078 and con-  
tained the active drug and Pharmatose® 110 M lactose monohydrate carrier  
20 untreated or polished. The mean particle size of the polished carrier was  
about 60 microns and there were no significant amounts of smaller than 40  
micron particles. The formulations were stored one week at 25°C/33 % RH  
and were then packed in two TAIFUN® MDPIs for testing the initial perform-  
ance of the products. Two polycarbonate tubes were filled with the same  
25 powder and placed immediately to stress circumstances of 45°C/75 % RH for  
one month. The tubes are permeable to moisture and do not shelter the  
formulation. Then two TAIFUN® MDPIs were filled with the formulation and  
tested.

30 The tests were done using an Andersen cascade impactor at constant ambi-  
ent circumstances 25 °C/60 % RH. The main parameter was fine particle

12

fraction, which is the percentage of smaller than 5.8  $\mu\text{m}$  drug particles of the total delivered dose. Each result is the mean of two tests. The dose strengths of the formulations were: salbutamol 50  $\mu\text{g}/\text{dose}$ , formoterol 12  $\mu\text{g}/\text{dose}$  and budesonide 100  $\mu\text{g}/\text{dose}$ . The results are seen in columns shown in figure 5. and explained below.

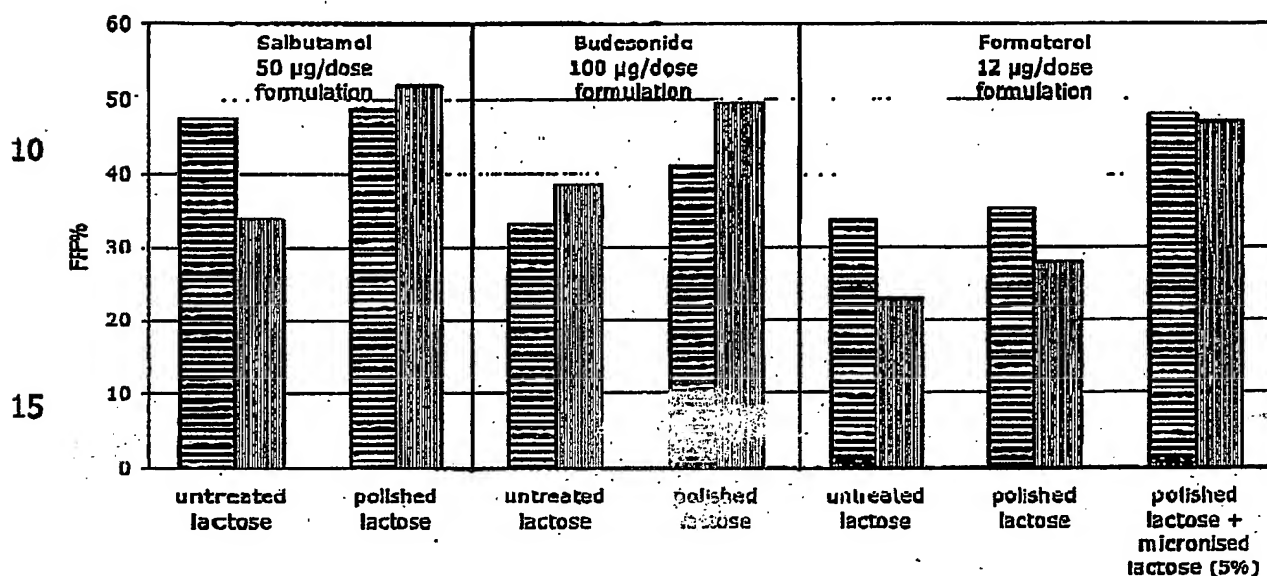

20 Figure 5. Fine Particle Fraction (FPF%)

#### 1. Salbutamol formulations

When untreated lactose was used, the initial FPF % is rather good (over 45 %), but considerable reduction to less than 35 % takes place during storage. For polished lactose, the initial value is almost 50 % and improves to over 50 % during storage.

## 2. Budesonide formulations

Initially the FPF value for untreated lactose is not good but improves somewhat during storage. For polished lactose the behaviour is very similar, but  
5 the FPF values are clearly better (average 35 contra 45 %). It is obvious that lipophilic budesonide stands stress better than hydrophilic salbutamol with some amorphous material in it. Also, no amorphous fraction was detected in budesonide in a microcalorimetric study. However, the effect of polished carrier is clear.

10

## 3. Formoterol

With untreated lactose the initial value is acceptable but falls down during storage. With polished lactose the result is clearly better, even if a small depletion in the FPF % is seen during stress. The best results were obtained,  
15 when a small amount of micronised carrier (5 w % of total) was mixed in the formulation. Now the FPF % was excellent and did not change when stressed.

20 It is commonly known that formoterol is difficult to formulate as inhalation powder; the FPF % is low and the physical stability is questionable. In this work it was detected that formoterol contains percents of amorphous material, which is difficult to recrystallise beforehand. By increasing the temperature for one hour to 60°C as hexane suspension (as described in patent FI  
25 105.078) most, but not all of the unstable material, recrystallised. This finding may explain the difficulties.

## POLISHING OPTIMIZATION

30 The polishing optimization parameters for lactose were determined with a device described earlier in the specification. The parameters involved were

rotation speed of the rotor, the distance between rotor and stator, amount of suspension medium (ethanol). the polishing time was 3 hours and the amount of lactose 400 g. The parameters used are listed in the table below:

## 5 POLISHING OPTIMIZATION

### PARAMETRES

	Test,	rpm (max 13000),	stator/rotor distance	amount of suspension medium
10	1.	60%	0.5mm	4000ml
	2.	25%	0.3mm	6000ml
	3.	60%	0.5mm	4000ml
	4.	60%	0.5mm	4000ml
	5.	25%	1.0mm	2000ml
15	6.	90%	0.3mm	2000ml
	7.	90%	0.3mm	6000ml
	8.	25%	1.0mm	6000ml
	9.	25%	0.3mm	2000ml
	10.	90%	1.0mm	2000ml
20	11.	90%	1.0mm	6000ml

Polishing time: 3 h; Amount of lactose: 400 g

The results obtained presented in Figure 6 reveal the rotation speed is the most important factor when high shape factors are desired.

25

Increase of micronsize material to the carrier is a well-known method to increase the FPF %. In EP 663.815 this method was utilised to increase FPF % when the formulation was used in a capsule based device. However, the flow properties of such a formulation are too impaired for use in multi-dose powder inhalers. Tests with Taifun MDPI showed that increase of micronised material to a formulation based on commercial lactose grades for Inhalation worsens the flow properties, followed by diminished dose and worse dose accuracy.

30

## 15

For acceptable dose accuracy, the total amount of micronised substance should not exceed 15 w-% of total. Polished carriers offer a new opportunity to benefit increase of micronised carrier to enhance both the FPF % and physical stability. Accordingly, the flow properties stay acceptable to be used in MDPIs, because the tested dose accuracy of the formulation with added micronsize carrier was excellent with an RSD % value of 7.

The method accepts different carrier materials, eg. glucose and mannitol were tested successfully in this work. Concerning the suspending liquid, the only prerequisite is that the solid material does not dissolve in the liquid. If the liquid is volatile (Bp less than 100°C), it can easily be dried in commercial dryers. If the liquid is less volatile, it can be washed with an appropriate volatile liquid during filtering. Total dryness is not needed if the carrier is used immediately in a wet mixing process for manufacture the final formulation, as explained in the Finnish patent FI 105 078. It is possible to have a certain fraction of polished particles in the final product. For example, a fraction of 20-30 µm polished particles can be returned to the main fraction. This can be done by selecting suitable filters to the filter pack so that the wanted fraction returns for retreatment. If no other means are used, the smallest particles, depending on the fine filter, will be present in the final product and are obviously located on the surface of the larger particles during the process or latest during drying. If the pore size of the finest filter is 2 microns, smaller carrier residues remain in the product and adhere on the surface of larger carrier particles during drying, preferably in a rotating evaporator. Also a soluble component may be added to the suspension for any reasons.

There are numerous ways to utilise polished carriers. As explained above, a combination of different fractions of polished carrier may be obtained for manufacture of the final formulation. The smaller size particles may act as ball bearings between the larger ones or they may form a separating layer on the surface of the larger particles. Then Van der Waals forces are dimin-



ished and disintegration of the drug particles is facilitated. The consequences can be easily studied by cascade impactor tests.

5 Polished particles may be coated with a secondary agent. Several methods are known including spray-drying with micronised or solvated secondary agent, gentle ball-milling with the agent (e.g. L-leucine and Mg-stearate) and gas-diffusion in vacuum. As reported, such coatings may greatly improve the flowability of carriers and also increase the FPF %. Polished carriers are excellent substrates for further developments.

10

The decisive idea in suspension-polishing is that vigorous modification of surfaces can be done without creating amorphous material. On the contrary, the surface layer with possible amorphous material is removed. The liquid acts as a coolant allowing only abrasion without melting or deforming the contact areas. The liquid also prevents agglomeration of particles of any size. If non-toxic liquids are used, toxicological hazards are avoided. A very important feature is that drastic enhancement in the properties can be done with physical treatment without using new chemical components, which should be proven safe for human inhalation before registration of the final medicinal product.

20

As the treatment can be carried out in a closed system using other than water-based liquids, microbial contamination can be avoided.

25 The method is ready to be scaled-up in manufacturing scale. The manufacturer of the used dispersers informs that the system can be enlarged to any scale. The principle of the method is so simple that no risks in scale-up are seen. Also the cross-flow filter is in industrial use in filtering waste water.

L 2

17

Claims

1. Method for treating a particulate carrier for an inhalation powder improving stability and flow properties of the carrier, **characterized** in that carrier  
5 is abraded suspended in a liquid medium into which the carrier is essentially insoluble, the liquid medium is removed and the carrier recovered.
2. Method according to claim 1, **characterized** in that the carrier is abraded  
10 with a mixing device using an effect below that required for crushing the carrier particles.
3. Method according to claim 1 or 2, **characterized** in that the rotation speed of the mixing device is lowered during the treatment.
- 15 4. Method according to any of claim 1 to 3, **characterized** in that the carrier suspension is recirculated to the mixer.
5. A method according to any of the proceeding claims, **characterized** in that the suspension is recirculated through a filter.  
20
6. A method according to claim 5, **characterized** in that a certain desired size range or ranges are recirculated to the mixing device.
7. A method according to any of the proceeding claims, **characterized** in  
25 that said media is a hydrocarbon, perfluorinated ether, fluorinated ether, perfluorinated hydrocarbon, fluorinated hydrocarbon, methanol, ethanol or any other alcohol or hydrocarbon.
8. A method according to any of the proceeding claims, **characterized** in  
30 that said carrier after filtration is used undried for formulation.
9. A method according to any of the proceeding claims, **characterized** in that said carrier is dried after filtration and stored for future used.

18

10. A method according to any of the proceeding claims, **characterized** in that the abraded carrier is at least partly covered particles smaller in size than said carrier.

5 11. A method according to claim 10, **characterized** in that the abraded carrier and the small sized particles are of the same material.

10 12. A method according to any of the proceeding claims, **characterized** in that the carrier to be abraded is lactose or a monohydrate thereof, glucose, mannitol, trehalose, sucrose, any other sugar, polysaccharide or any other compound used as a carrier.

15 13. Carrier for an inhalation powder, which carrier is stable and possesses good flowing properties, **characterized** in that the carrier is abraded suspended in a liquid medium, in which said carrier is essentially insoluble.

20 14. Carrier according to claim 13, **characterized** in that that the carrier is abraded with a mixing device using an effect below that required for crushing the carrier particles.

25 15. Carrier according to claim 13 or 14, **characterized** in that the carrier is filtrated and used for formulation undried or dried and stored for future use.

30 16. Carrier according to any of the claims 13 - 15, **characterised** in that the filtrated carrier contains more than one main range of particle sizes of abraded carrier.

17. Carrier according to any of the proceeding claims, **characterized** in that the carrier to be abraded is lactose or a monohydrate thereof, glucose, mannitol, trehalose, sucrose, any other sugar, polysaccharide or any other compound used as a carrier.


## 19

18. Preparation for Inhalation purposes comprising an active agent, a carrier and optional excipients used in inhalable preparation, **characterized** in that at least a part of the carrier used is abraded suspended in a liquid medium, in which the carrier is essentially insoluble.

5

19. A preparation according to claim 18, **characterized** in that carrier contains more than one main range of particle sizes.

Jos lisätään asiaa 14/15 väliin, pitäisikö ko. asia mainita myös tässä ?

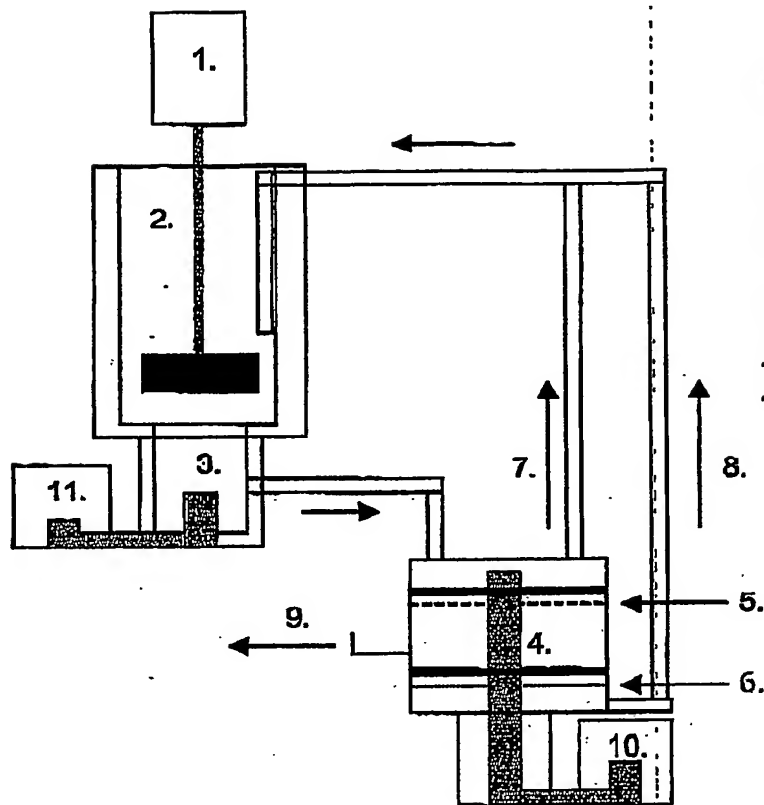

1. Mechanical mixer
2. Water cooled suspension vessel
3. Super Dispax™ water cooled stator/rotor chamber
4. Cross flow classifying filter pack
5. Coarse filter ( 40 microns )
6. Fine filter ( 0.5 microns )
7. Return, coarse particles
8. Return, supernatant
9. Output ( particles between filters )
10. Filter motor
11. Super Dispax™ motor

FIG. 1

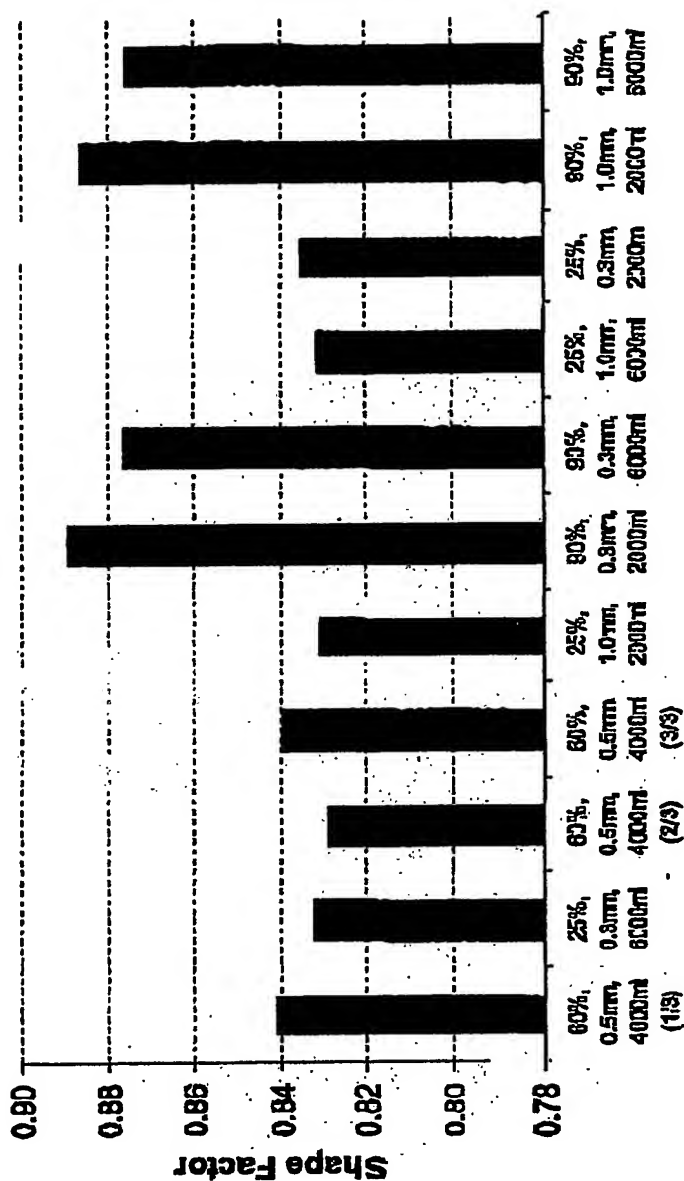


FIG. 6

**This Page is Inserted by IFW Indexing and Scanning  
Operations and is not part of the Official Record**

**BEST AVAILABLE IMAGES**

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

☒ **BLACK BORDERS**

☐ **IMAGE CUT OFF AT TOP, BOTTOM OR SIDES**

☐ **FADED TEXT OR DRAWING**

☒ **BLURRED OR ILLEGIBLE TEXT OR DRAWING**

☐ **SKEWED/SLANTED IMAGES**

☐ **COLOR OR BLACK AND WHITE PHOTOGRAPHS**

☐ **GRAY SCALE DOCUMENTS**

☒ **LINES OR MARKS ON ORIGINAL DOCUMENT**

☐ **REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY**

☐ **OTHER:** \_\_\_\_\_

**IMAGES ARE BEST AVAILABLE COPY.**

**As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.**